



PRIMEKINETIX™

INNOVATIVE MOVES

MyoStim Maxx

USER MANUAL

Contents

Cautions-----	01
Attachment list-----	02
Functional area -----	03
How to stick the gel pad-----	03
Charging Instructions-----	04
Using-----	04
Massage Host Operation Instructions-----	05
Maintenance instruction-----	06
Storage-----	06
Specifications-----	07
Simple troubleshooting-----	07
Description of hazardous substances-----	08

Cautions

- Please use product after reading the user manual carefully, and keep the manual properly in case you need it at any time.
- Please keep away from children.
- Product should be attached to the skin first before using, and then turn on.
- When using, please take off electrical conductors to avoid damage or Static electricity, such as the electronic watch, finger rings, earrings.
- Please stop using the product if you feel unwell during use.
- Excessive use of the same part may cause skin redness, please stop using.
- Please be sure to use our brand gel pad provided by us, and recommend each set of gel pad only for one user.
- The disposal of unit and battery should be classified.




Compulsory Compliance

- Do not prick or tear the product hard.
- Keep away from water and wet environment in case of problem.
- Do not take off the massage pad or change the massage position when the power is on.
- Do not use on head, heart area, tattoo area or the genital area.
- Do not use device while driving, strenuous exercising, sleeping or bathing.
- Please do not use on children or language disorders.
- The elderly and those who cannot take care of themselves need to be accompanied when using.
- Do not store in direct sunlight, dusty places, high temperature or humidity.
- Do not disassemble, install or repair.
- Do not use device with other medical devices at same time.


Forbidden and prudent people

- Do not use after excessive drinking.
- Do not use immediately on the abdomen after meals.
- Acute patients, blood patients and people with severe arrhythmia, cardiovascular and cerebrovascular diseases, malignant tumor, epilepsy, disturbance of perception.
- Do not use on scytitis, skin sensitivity, inflammation of the breast, burns, bleeding, itching area.
- Do not use in pregnancy, during menstruation, children under 12 years old.
- People who are allergic to the gels are not allowed to use this device.
- People who implanted electronic items in the body.(such as heart pacemaker, heart rate regulator, heart stent).
- On some kind of therapy or symptoms began with fever, make sure to consult your doctor before using.
- Patients who have fever or are receiving treatment, please follow the doctor's advice.
- Patients with severe abnormal blood pressure, severe anemia and hypoglycemia, please follow the doctor's advice.


Attachment list

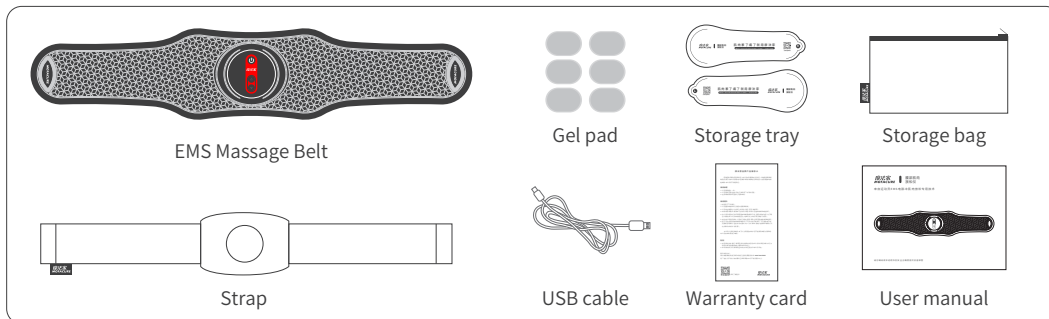


Functional area

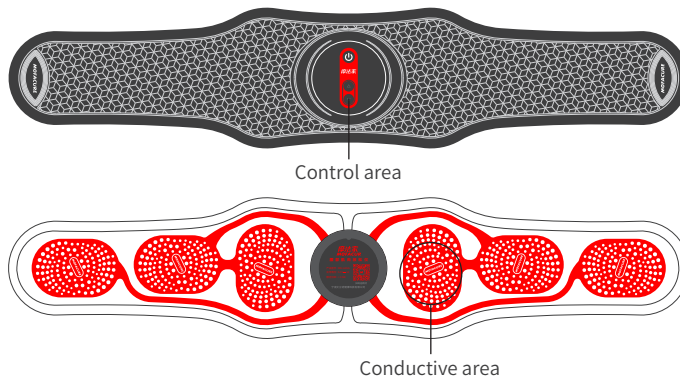


How to stick the gel pad

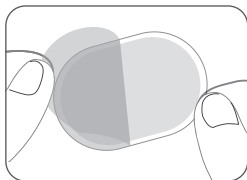

1

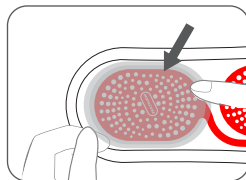


2

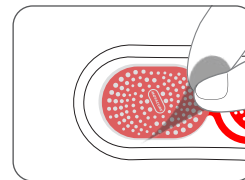

3

1. Take out the conductive gel pad, remove the blue plastic film from one side of the gel pad.
 2. Follow the schematic, stick the gel pad to the conductive area.
 3. Remove the transparent film of the gel pad, then stick the product to the body and turn it on.
- (note: Please try to avoid the skew, wrinkle and bubble of gel pad. If there are bubbles when sticking the gel pad, please gently push the bubbles to the transparent plastic film.)

Charging Instructions

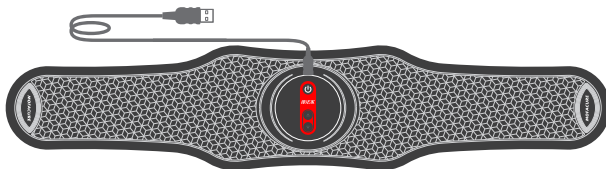

Before use, it is recommended to use the attached USB cable to charge the device. The USB cable can be connected to the computer or adapter, and the charging time is 2 hours.

State of indicator light during charging: when charging, the red indicator light flashes; when fully charged, the red indicator light is always on.

Standard charging voltage and charging current: 5.0V / 1.0A


(note: To avoid damaging internal devices, it is forbidden to use adapter with charging voltage greater than 5V and charging current greater than 1A to charge the equipment)

Using

Before using, please keep the body parts clean and dry, then peel off the transparent film of the gel pad. Place the device onto the treatment areas(Waist, Belly). Please press down firmly to ensure a full and firm contact with the skin. Please start from the low intensity at the first time to use, and then gradually increase the intensity.



Abdomen

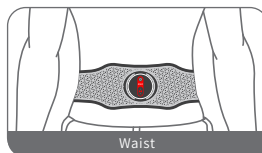

Waist

Three modes Five massage techniques Ten combination frequency-programming

(Use Mofacure for muscle ache and pain)



Beating



Knocking



Kneading



Manipulation



Massage






15 intensity levels



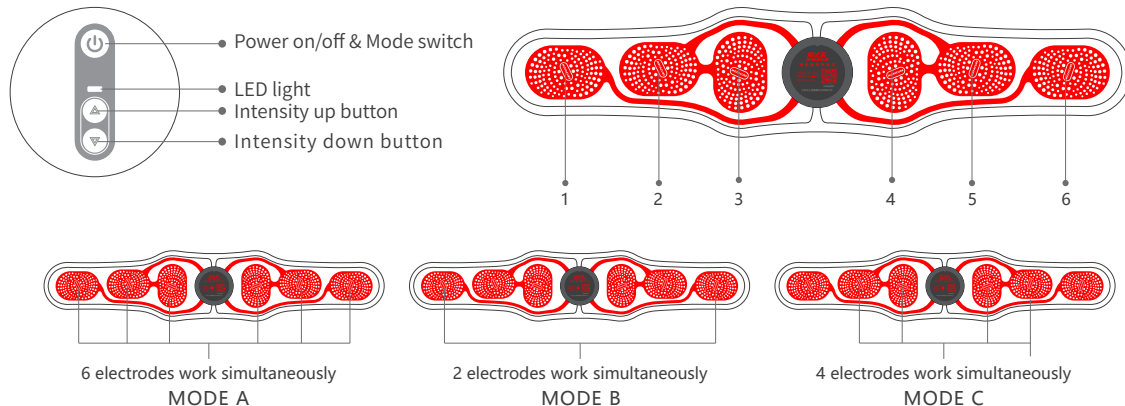

20minutes timer

Massage Host Operation Instructions



Long Press the  button to turn on the device, the green indicator flashes once and it beeps once. Long press the  button to turn off the device, the it beeps once.

Press the "" or "" button to adjust the intensity, the output intensity total have 15 levels, when turn on the device, the intensity is on level 1 .When you adjust the intensity the green indicator flashes once and it beeps once .When it's in the strongest intensity or lowest intensity,the green indicator flashes twice and it beeps twice as a reminder.

When turn on the device, the program defaults to mode A with level 1, electrode 1,2,3,4,5,6 will work at the same time. If press the power button again, it will be changed to mode B and electrode 1, 6 will work at same time. Press the power button again, it will be changed to mode C and electrode 2,3,4,5 will work at same time. The three modes will work in cycle operation.


When the device is working, if it be removed from the skin, or the device is no stick close to skin tightly,the green indicator continuously flashes and it continuously beeps for one minute as a reminder, and then automatically turn off. When the device return to stick the skin again, the intensity will return to level 1.

Maintenance instruction

- Before using it at the first time, please charge the battery first. To extend the battery life, please do not use fast charging adapters. Get it fully charged before long-term storage and charge it every two or three months to extend the life of built-in battery.
- Before use, wipe the part that you want to massage with wet towel. and remove oil stain and sweat stain, to extend the service life of gel pad.
- After use, please gently remove one end of the device, to avoid gel pad sticking together, extend the service life of gel pad.
- Generally, gel pad can be used repeatedly for about 30-50 times. If the viscosity of the gel is reduced, dip in water with fingers and then stroked the gel pad surface to make it moist. After drying, it will restore some viscosity.
- If gel pad is not sticky enough, it cannot reach the effect of massage, and it is easy to cause tingling sensation. If the adhesive properties of the gel cannot be restored, please replace a new one.
- If there is a stain on the surface of the product, please gently wipe it with a wet towel. If the gel has stains, dip in water with fingers and then stroked the gel pad surface to make it moist and clean.
- Please keep product away from the following environment: in the sun, in the place of strong electromagnetic field, hot and humid, dusty, electromagnetic environment.

Storage

- As picture: after remove from the body, please stick it to the storage board for storage.
- Please do not let the gels touch each other, such as folding the product after the storage tray is not pasted.
- Please do not fold the storage board, otherwise it will cause deformation and cannot be used normally.
- If the waist protector is to be put into the storage bag, please put the sponge in the middle of the belt, and then bend the belt and insert it into the storage bag to prevent deformation and indentation caused by extrusion.

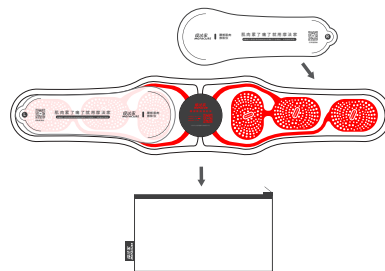

Safety caution:

- Please keep the storage tray clean;
- Do not place the product in high temperature, damp places, avoid direct sunlight, (such as in the car).
- Please keep away from children and pets, place device in a safe place.

Specifications

Name	Myostim Maxx	Weight	108g	Frequency	1-100Hz
Material	Silicone、PC+ABS	Pulse width	200μs	Operating temperature	0°C-45°C
Model No.	MFJ-6302L	Output Voltage	152 V at 1K Ω	Operating humidity	Relative humiditys ≤85%
Size	580x120x16mm	Output power	< 2.5 mW	Battery Capacity	185mAh Rechargeable Polymer Lithium battery

Simple troubleshooting

Device cannot turn on	The battery is low, please charge the product.
Abnormal shutdown	The product does not stick close to the skin, please stick to the skin tightly or wet the place with towel before turn it on.
Massage intensity is too weak	1. If low intensity, please press the key of "+". 2. If poor viscosity of gel, please replace a new one.
Cannot be charged	Please make sure the adapter Output Voltage is DC 5.0V / 1.0A and in normal working condition.
Sting feeling during massage	1.Please ensure that the gel pad stick smoothly with the skin. If the gel is not sticky enough, clean the wet gel with clean water, use it after restoring the viscosity, or replace the new gel pad. 2.If you feel tingling at any position during each use, you may be sensitive, and this product is not applicable.

Description of hazardous substances

	Shell	PCB	Battery	Electrode Slice	USB Cable	Package
Pb	○	○	○	○	○	○
Hg	○	○	○	○	○	○
Cd	○	○	○	○	○	○
Cr, VI	○	○	○	○	○	○
PBB	○	○	○	○	○	○
PBDE	○	○	○	○	○	○

O: Indicates that this toxic or hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement in GB/T 26572-2011.

X: Indicates that this toxic or hazardous substance contained in at least one of the homogeneous materials for this component is above the limit requirement in GB/T 26572-2011.

All components with "X" in the table are in accordance with the EU ROHS recast regulation 65 / 2011 - / EU directive.